



HEAT SEALED POWDER COAT STEEL PASS THROUGH 20 SCOPE CAPACITY CABINET, KEY LOCK

SC80WTD RDP-20

*Also available
in a Stainless Steel Finish*

FEATURES

Scopes not
included

- **Double sided pass through cabinet allows access between both procedure and reprocessing rooms**
- In-wall system with multiple available locking options
- Positive pressure fully redundant fan system
- Industries only unit with 4 intake and 2 exhaust fans
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Easy to clean removable dual drip trays
- Holds 20 scopes with our secure adjustable mounting system
- Replaceable HEPA Filter system without removal of scopes
- Medical grade electrical system
- Tempered glass doors, one on each side
- Cabinets delivered fully assembled
- Cabinet dimensions: 93.25" H x 41" W x 28" D
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada



COLORS

Available in more than twenty color options.

View them online: <https://www.harloff.com/finish-options/>

Simple "Plug & Play" upgradeable with Dri-Scope Aid®, industry leading internal lumen drying system

HARLOFF
INNOVATING STORAGE FOR CARE

Engineered &
Manufactured
in the U.S.



View similar products online at
<https://www.harloff.com/suredry-scope-cabinets/>

**Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300**

The SureDry Scope Cabinet includes the industry best Positive Pressure HEPA Filtered Drying system. Inside air is replaced inside the cabinet four times per minute.


**Optional dimmable LED
lighting also available!**

Why is a positive pressure system important?

The unique design of the SureDry Scope Cabinet includes both top fans for intake and lower fans for exhaust. This push/pull design creates a vortex of air that flows through the cabinet ensuring maximum drying efficiency.

Intake: Intake airflow comes through separate HEPA filters (which remove 99.7% of particulates, .03 Microns) driven by fully redundant fans producing a positive pressure vortex from the cabinet top. This creates an even circular motion of airflow over all scopes.

Exhaust: The lower exhaust fans pull the circulated air out of the cabinet ensuring the flow of air is consistent from top to bottom, providing a complete drying system for the entire scope. This push/pull multiple fan system (an industry first) ensures a positive pressure system to meet all major guidelines and it ensures forced drying and equal drying levels across the full length of each and every scope.

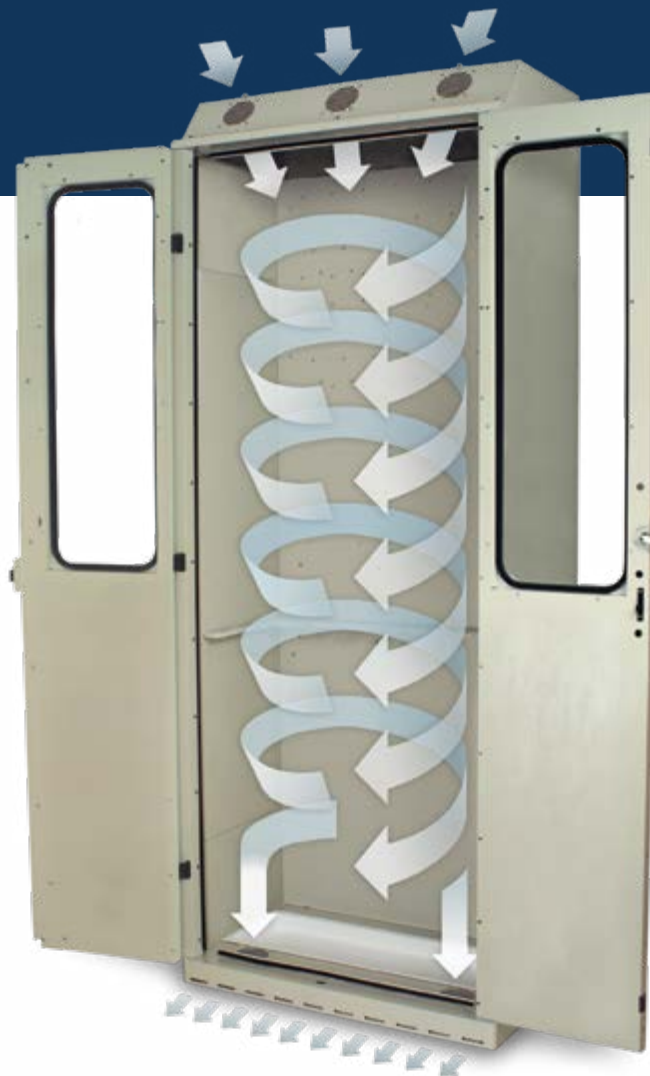

ST91:2021

SureDry Cabinets are a Smart Choice to Meet and Exceed ANSI/AAMI ST91:2021 and Other Guidelines

The references below are taken from these guidelines. Channel drying available only with the addition of the Dri-Scope Aid® system.

☒ **Scopes Stored Vertically**

11.2.1 "stored vertically, the insertion tube should be as straight as possible, with the distal tip hanging freely."
11.2.1 "When flexible and semi-rigid endoscopes are hung in the vertical position, coiling or kinking is prevented."

☒ **Hanging Freely Without Coiling**

11.2.2 "should be stored in a cabinet that is of sufficient height, width, and depth to allow flexible endoscopes to hang vertically without coiling and without any components touching the bottom of the cabinet."

☒ **Channel Drying + Forced Air**

8.2.5.1 "Flexible endoscopes with channels should be dried for a minimum of 10-minutes with pressure-regulated forced instrument air or a minimum of HEPA-filtered air." | 8.2.5.2 Use pressure-regulated instrument air or HEPA-filtered air to dry the channels in accordance with the manufacturer's written IFU.

☒ **HEPA Filtered Positive Pressure**

11.2.2.2 Cabinets that enable circulation of HEPA-filtered or instrument air through the cabinet at... positive pressure.
NOTE: SureDry cabinets feature an industry leading 3 intake and 2 exhaust fan system.

☒ **Easy Scope Identification***

11.2.3 "Protocols should be developed to help ensure that users can readily identify an endoscope that has been processed and is ready for patient use." * Available with Tracker-equipped SureDry scope cabinets

☒ **Single-Use Connectors**

Use of single use connectors ensures no cross-contamination as well as simplifies staff efficiency.

☒ **Clinical Documentation**

ASGE "Drying the interior of an endoscope should use a sufficiently prolonged flow of medical air through all accessible channels, ideally simultaneously and for at least 10 minutes, for greatest efficiency." Per the study, "Multisociety guideline on reprocessing flexible GI endoscopes: 2016 update."
NOTE: SureDry cabinets are the only clinically proven vertical drying cabinet.