

DSC36T GK-DP *with Glass Windowed Door*

TALL 36 DILATOR HEPA DRYING CART for Savary-Gilliard® Dilators up to 39" in length


Current guidelines require both HEPA filtration and that endoscopes be hung vertically due to the need for inner channels to drain all excess fluid after cleaning. Multiple facilities have requested a cart or cabinet to allow them to store their Savary-Gilliard® dilators vertically because they also have an internal channel.

With this in mind, the design and engineering team at Harloff has introduced the new Endoscopy Dilator Storage Cart. Using the industry leading design of our Scope Drying cabinets we produced a HEPA filtered Positive Pressure drying system with a fully adjustable fan that provides Savary dilator vertical drying for up to 36 dilators 39" or less in length.

Features:

- Positive pressure fan system dries the full Savary-Gilliard® dilators and other wire guided dilators
- Includes six Savary pull out slides that each hold six dilators for easy access to similar sized dilators at one time
- Replaceable HEPA Filter removes 99.7% of particulates, .03 Microns
- Fan run time (from 10 to 180 minutes) is fully adjustable and activated each time the dilator storage cart is closed
- Durable scratch resistant heat sealed powder coat finish
- Each dilator storage cart can hold up to 36 individual dilators, varying in size from 5mm (15 French) to 20mm (60 French) and up to 39" in length
- Choice of lock options: keyless, latch lock, standard key lock or standard electronic lock
- Medical grade electrical system
- Sturdy, fully welded steel body
- Industry best 5" quiet rolling casters - one steering caster and two braking casters
- Industry only fully curved protective bumper

Color shown is
Hammer-tone Gray.
More than 20 additional
colors are available.

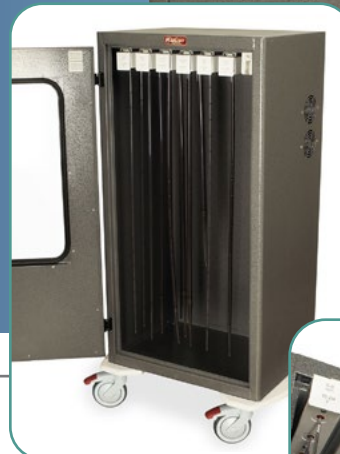

Engineered &
Manufactured
in the U.S.

Assembled Cart Dimensions:

- 51.88" H x 29.88" W x 22" D (131.8 cm H x 75.9 cm W x 55.9 cm D)
- Cart weight: 109 lbs (49.5 kg)
- Approximate shipping weight: 134 lbs (60.8 kg)
- Approximate Shipping Dims: 50" H x 34" W x 24" D (127 cm x 86.4 cm x 61 cm)

Warranty: 12-Year Manufacturer's Warranty on shipments to US and Canada, 5-Year Manufacturer's Warranty on all other International shipments. 1-Year Limited Warranty on any bought-out parts

HARLOFF
INNOVATING STORAGE FOR CARE

View similar products online at
<https://www.harloff.com/product-category/dilator-carts/>

Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300