



## STAINLESS STEEL 10 SCOPE CAPACITY CABINET, KEY LOCK

Scopes not included

**SCSS8030DRDP-DSS2310**

### FEATURES

- Positive pressure fully redundant fan system
- Unit features two intake fans and one exhaust fan
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter system removes 99.7% of particulates, .03 Microns – filters are easily replaced without requiring scope removal
- Scope capacities of 10, 8, 5 or a custom setup using our secure adjustable mounting system
- Configurations available for bronchoscopes, cystoscopes, TEE Probes and many other scopes
- Medical grade electrical system
- Tempered glass double doors
- Multiple locking options, including badge entry
- Easy to clean removable drip tray
- Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada
- Cabinet dimensions: 93" H x 30" W x 24" D
- Stainless steel is commonly used in the manufacturing of medical devices due to its rust and corrosion resistance, anti-bacterial properties and ease of cleaning and sterilization.



**Includes Dri-Scope Aid®, industry leading internal lumen drying system**

Please refer to Dri-Scope Aid® specification sheet for more information



Engineered &  
Manufactured  
in the U.S.



View similar products online at  
<https://www.harloff.com/suredry-scope-cabinets/>

**Toll Free: (800) 433-4064 • [medicalsales@harloff.com](mailto:medicalsales@harloff.com)  
[www.harloff.com](http://www.harloff.com) • Ph: (719) 637-0300**

The ultimate in Endoscope drying systems is now available! SureDry Scope Cabinets feature an external HEPA filtered positive pressure drying system. Dri-Scope Aid® is the industry leading internal lumen drying system. This combination provides facilities with the assurance that they will be able to meet all current and anticipated future guidelines for external and internal drying of endoscopes, providing their patients the best possible assurance of infection prevention.

## DRI-SCOPE INTERNAL LUMEN DRYING SYSTEM



### FEATURES & BENEFITS

- Significantly improves drying of internal channels
- HEPA filter offers a 99.97% retention rate of particles that have a size > 0.3µm
- Provides constant air delivery through all internal endoscope channels
- Short drying cycle for quick turn-around time
- Continuous venting cycle for overnight storage
- Automatic high-pressure shutoff
- Fulfills forced air recommendations
- Self-adjusting high-pressure regulation

### EASE OF OPERATION

- Simple to attach individual scope connectors
- Hands-free operation with visual display of purging times
- Works with all major endoscope brands
- Disposable connectors to meet SGNA and AORN guidelines
- Dries 1 to 10 endoscopes

### OVERALL CABINET DIMENSIONS

- 93" H x 30" W x 24" D



- Constant Flow
- Consistent Pressure
- Repeatable Results

The Dri-Scope Aid® Cabinet supplies HEPA-filtered air to all internal channels to aid in drying endoscopes with an automated drying cycle for improved turn-around time, and a continuous venting cycle for overnight storage. The Dri-Scope Aid® system in the SureDry Scope Cabinets provides a complete drying system for your facility.

#### OPTIONAL: Customizable Hang Time Monitor

Adjust the number of days your scopes hang. You will be alerted when a scope is about to expire



Engineered &  
Manufactured  
in the U.S.



View similar products online at  
<https://www.harloff.com/suredry-scope-cabinets/>

Toll Free: (800) 433-4064 • [medicalsales@harloff.com](mailto:medicalsales@harloff.com)  
[www.harloff.com](http://www.harloff.com) • Ph: (719) 637-0300

**Optional dimmable LED  
lighting also available!**

# Why is a positive pressure system important?

The unique design of the SureDry Scope Cabinet includes both top fans for intake and lower fans for exhaust. This push/pull design creates a vortex of air that flows through the cabinet ensuring maximum drying efficiency.

## Intake:

Intake airflow comes through separate HEPA filters (which remove 99.7% of particulates, .03 Microns) driven by fully redundant fans producing a positive pressure vortex from the cabinet top. This creates an even circular motion of airflow over all scopes.

## Exhaust:

The lower exhaust fans pull the circulated air out of the cabinet ensuring the flow of air is consistent from top to bottom, providing a complete drying system for the entire scope. This push/pull multiple fan system (an industry first) ensures a positive pressure system to meet all major guidelines and it ensures forced drying and equal drying levels across the full length of each and every scope.



## SureDry Cabinets: A Smart Choice

to Meet and Exceed Current AORN, SGNA and AAMI Guidelines

Scope Storage Guidelines	Forced Air Drying	HEPA Filtered Air	Hang Scopes Vertically	Scopes Must Hang Freely, No Touching or Coiling	Easy to Clean and Disinfect	Closed Door Cabinet	Positive Air Pressure	Circulated Air Environment	Connectors — Use Disposables or Reusables Stored w/Scopes	Dry Scopes Completely	Uses Channel Drying (with Dri-Scope Aid®)
ANSI/AAMI ST91:2021		SureDry scope cabinets, with or without Dri-Scope Aid®, meet and exceed updated guidelines									
AORN	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
SGNA	✓	✓	✓	✓	✓	no data	✓	no data	✓	✓	✓
AAMI	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓



# 10 Reasons to Choose **HARLOFF** for Your Scope Cabinets

*A Smart Choice to Meet & Exceed Current  
SGNA, AORN and AAMI Guidelines*



**1** *Industry Best Positive Pressure HEPA  
Filtered Drying System.* Designed to add Dri-Scope  
Aid® now or later with easy “Plug & Play” installation.

## **Proudly Made in the USA**



**6** Harloff products are proudly  
*Designed, Engineered and  
Manufactured in the United States*  
at our Colorado Springs, CO facility.




## **LED Lighting Option**

**2** *Brighten your cabinet*  
and more easily inspect  
scopes with fully dimmable  
optional overhead LED lighting  
which triggers on/off when you  
open and close the door.

## **Easy Access to HEPA Filters**

Change HEPA filters  
*without opening the  
cabinet.* Save time and  
ensure no cycle interruptions.

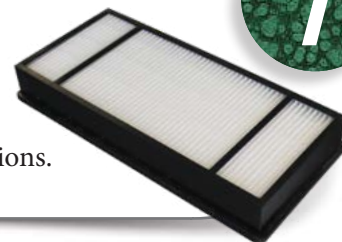

## **Disposable Connectors**

*Save time and ensure  
compliance!* Use single-use  
connectors for all internal  
channels of the endoscope.  
Part of the Dri-Scope Aid® package.



## **Keyed or Keyless Locks**

*Options include* key locks,  
basic electronic push  
button locks, advanced  
keypads and proximity  
readers. View locking options online  
at [www.harloff.com/locking-options/](http://www.harloff.com/locking-options/)



## **Hooks and Holders for Specialized Scopes**

Major brands to include  
*Fuji, Olympus, Pentax, Karl Storz and others.*



## **Finish Options and Colors**

Stainless Steel, Heat Sealed Powder Coat  
(High Gloss or Textured Low Gloss inside & out,  
view online at [www.harloff.com/finish-options/](http://www.harloff.com/finish-options/))  
or upgrade to an *antimicrobial finish.*

## **12 Year Warranty**

*Our Manufacturer's  
Warranty is Unmatched* in the  
industry. We build our cabinets to  
last and we stand behind them.

## **Mobile Scope Cabinet Kit**

Convert any size of scope cabinet  
into a mobile version. Includes 4  
3" casters and steel mounting plate.

