### STAINLESS STEEL 15 SCOPE CAPACITY CABINET, KEY LOCK



SCSS8036DRDP-DSS3316

Scopes not included

### **FEATURES**

- · Positive pressure fully redundant fan system
- Industries only unit with 3 intake and 2 exhaust fans
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Holds 15 scopes with secure adjustable mounting system
- Replaceable HEPA Filter system without removal of scopes
- · Medical grade electrical system
- · Tempered glass double doors
- Multiple locking options available including badge entry
- · Easy to clean removable drip tray
- · Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada
- Overall cabinet dimensions: 93"H x 36"W x 21"D
- Stainless steel is commonly used in the manufacturing of medical devices due to its rust and corrosion resistance, anti-bacterial properties and ease of cleaning and sterilization.



Includes Dri-Scope Aid®, industry leading internal lumen drying system

Please refer to Dri-Scope Aid® specification sheet for more information



Engineered & Manufactured in the U.S.



View similar products online at

https://www.harloff.com/scope-drying-and-storage-cabinets/suredry-scope-cabinets/

Toll Free: (800) 433-4064 • medicalsales@harloff.com www.harloff.com • Ph: (719) 637-0300 The ultimate in Endoscope drying systems is now available! SureDry Scope Cabinets feature an external HEPA filtered positive pressure drying system. Dri-Scope Aid® is the industry leading internal lumen drying system. This combination provides facilities with the assurance that they will be able to meet all current and anticipated future guidelines for external and internal drying of endoscopes, providing their patients the best possible assurance of infection prevention.

### **DRI-SCOPE INTERNAL LUMEN DRYING SYSTEM**



- Constant Flow
- Consistent Pressure
  - Repeatable Results

The Dri-Scope Aid® Cabinet supplies

automated drying cycle for improved

Dri-Scope Aid® system in the SureDry

Scope Cabinets provides a complete

drying system for your facility.

turn-around time, and a continuous venting cycle for overnight storage. The

HEPA-filtered air to all internal channels to aid in drying endoscopes with an

### **FEATURES & BENEFITS**


- Significantly improves drying of internal channels
- HEPA filter offers a 99.97% retention rate of particles that have a size > 0.3µm
- Provides constant air delivery through all internal endoscope channels
- · Short drying cycle for quick turn-around time
- Continuous venting cycle for overnight storage
- · Automatic high-pressure shutoff
- · Fulfills forced air recommendations
- · Self-adjusting high-pressure regulation

#### EASE OF OPERATION

- · Hands-free operation
- · Visual display of purging times
- Simple to attach individual scope connectors
- · Works with all major endoscope brands

INNOVATING STORAGE FOR CARE

· Filter replacement indicator light

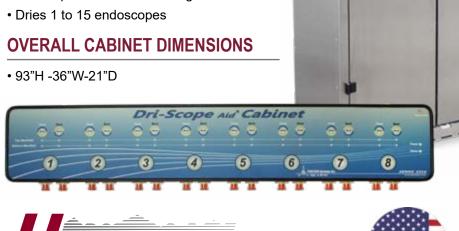

Engineered & Manufactured in the U.S.



View similar products online at

https://www.harloff.com/scope-drying-and-storage-cabinets/suredry-scope-cabinets/

Toll Free: (800) 433-4064 • medicalsales@harloff.com www.harloff.com • Ph: (719) 637-0300 Why is a positive pressure system important?

The unique design of the SureDry Scope Cabinet includes both top fans for intake and lower fans for exhaust. This push/pull design creates a vortex of air that flows through the cabinet ensuring maximum drying efficiency.

#### Intake:

Intake Airflow comes through 3 Separate HEPA filters (removes 99.7% of particulates, .03 Microns) driven by 3 fully redundant fans producing a positive pressure vortex from the top of the cabinet. This creates a circular motion of airflow evenly over all scopes.

### **Exhaust:**

The lower exhaust fans pull the circulated air out of the cabinet ensuring the flow of air is consistent from top to bottom, providing a complete drying system for the entire scope. This push pull multiple fan system, an industry first, ensures not only a positive pressure system to meet all major guidelines but also ensures forced drying and equal levels across each and every scope across it's full length.



A Sm	A Smart Choice to Meet & Exceed Current AORN, SGNA and AAMI Guidelines												
Scope Storage		Forced Air Drying	HEPA Filtered Air	Hang Scopes Vertically	Scopes Must Hang Freely, No Touching or Coiling	Easy to Clean and Disinfect	Closed Door Cabinet	Positive Air Pressure	Circulated Air Environment	Connectors — Use Disposables or Reuse- ables Stored w/Scopes	Dry Scopes Completely	Use Channel Drying*	
AOR	2N	no data	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	no data	no data	
SGN	Α	$\checkmark$	$\checkmark$	<b>✓</b>	<b>✓</b>	<b>✓</b>	no data	no data	no data	<b>✓</b>	$\checkmark$	no data	
AAI	ΛI	<b>✓</b>	<b>√</b>	<b>√</b>	<b>✓</b>	<b>✓</b>	<b>✓</b>	<b>✓</b>	<b>✓</b>	<b>✓</b>	<b>√</b>	$\checkmark$	

Harloff Scope Cabinets Meet & Exceed These Guidelines (\*Channel Drying Available Only with Dri-Scope Aid®)

# 10 Reasons to Choose **HARLOFF** for Your Scope Cabinets





**Industry Best Positive Pressure HEPA** Filtered Drying System. Designed to add Dri-Scope

Aid® now or later with easy "Plug & Play installation.

A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines

### Proudly Made in the USA




Harloff products are proudly Designed, Engineered and Manufactured in the United States at our Colorado Springs, CO facility.


## Dri-Scope Aid

### Internal Lumen Drying System

HEPA filtered air to all internal channels on DSS3316 models only. Automated drying & continuous venting cycle for overnight storage. All SureDry cabinets are fully retrofittable with Dri-Scope Aid®.

## Easy Access to HEPA Filters

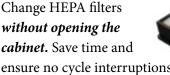

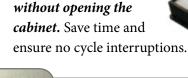

### Disposable Connectors

Save time and ensure compliance! Use single-use connectors for all internal channels of the endoscope. Part of the Dri-Scope Aid® package.

## **Keyed or Keyless Locks**

Options include key locks, basic electronic push button locks, advanced keypads and proximity

readers. View locking options online at www.harloff.com/locking-options/



## Hooks and Holders for **Specialized Scopes**

Major brands to include Fuji, Olympus, Pentax, Karl Storz and others.

## **Scope Cabinet Finish Options**

Choose Stainless Steel, Heat Sealed Powder Coat in High Gloss or Textured Low Gloss, or upgrade to an antimicrobial finish.

### **12 Year Warranty**

Our Manufacturer's Warranty is Unmatched in the industry. We build our cabinets to last and we stand behind them.



### **Superior Color Choices**

Over 20 Colors to choose from. Cabinets have a heat sealed powder coat online at <a href="https://www.harloff.com/finish-options/">www.harloff.com/finish-options/</a>