

# STAINLESS STEEL HIGH VOLUME 16 SCOPE CAPACITY CABINET, KEY LOCK



**SCSS8044DRDP-DSS3316**

Scopes not  
included

Engineered &  
Manufactured  
in the U.S.



## FEATURES

- *Industries most spacious scope cabinet allows for easy hanging and retrieval of all scope types*
- Positive pressure fully redundant fan system
- Industries only unit with 3 intake & 2 exhaust fans
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter system removes 99.7% of particulates, .03 Microns – filters are easily replaced without requiring scope removal
- Scope capacities of 16, 8 or a custom setup using our secure adjustable mounting system
- Configurations available for bronchoscopes, cystoscopes, TEE Probes & many other scopes
- Medical grade electrical system
- Tempered glass double doors
- Multiple locking options available including badge entry
- Easy to clean removable drip tray
- Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada
- Cabinet dimensions: 93" H x 44" W x 24" D
- Stainless steel is commonly used in the manufacturing of medical devices due to its rust and corrosion resistance, anti-bacterial properties and ease of cleaning and sterilization.



**Includes Dri-Scope Aid®, the  
industry leading internal lumen  
drying system**

The ultimate in Endoscope drying systems is now available! SureDry Scope Cabinets feature an external HEPA filtered positive pressure drying system. Dri-Scope Aid® is the industry leading internal lumen drying system. This combination provides facilities with the assurance that they will be able to meet all current and anticipated future guidelines for external and internal drying of endoscopes, providing their patients the best possible assurance of infection prevention.

View similar products online at  
<https://www.harloff.com/suredry-scope-cabinets/>

**HARLOFF**  
INNOVATING STORAGE FOR CARE

**Toll Free: (800) 433-4064**  
**medicalsales@harloff.com**  
**www.harloff.com**  
**Ph: (719) 637-0300**



## DRI-SCOPE INTERNAL LUMEN DRYING SYSTEM

- Constant Flow
- Consistent Pressure
- Repeatable Results

### FEATURES & BENEFITS

- Significantly improves drying of internal channels
- HEPA filter offers a 99.97% retention rate of particles that have a size  $> 0.3\mu\text{m}$
- Provides constant air delivery through all internal endoscope channels
- Short drying cycle for quick turn-around time
- Continuous venting cycle for overnight storage
- Automatic high-pressure shutoff
- Fulfills forced air recommendations
- Self-adjusting high-pressure regulation

### EASE OF OPERATION

- Simple to attach individual scope connectors
- Hands-free operation with visual display of purging times
- Works with all major endoscope brands
- Disposable connectors to meet SGNA and AORN guidelines
- Dries 1 to 16 endoscopes

### OVERALL CABINET DIMENSIONS

- 93" H x 44" W x 24" D

The Dri-Scope Aid® Cabinet supplies HEPA-filtered air to all internal channels to aid in drying endoscopes with an automated drying cycle for improved turn-around time, and a continuous venting cycle for overnight storage. The Dri-Scope Aid® system in the SureDry Scope Cabinets provides a complete drying system for your facility.



**OPTIONAL:**  
**Customizable Hang Time Monitor**  
Adjust the number of days your scopes hang. You will be alerted when a scope is about to expire



Engineered & Manufactured in the U.S.



View similar products online at  
<https://www.harloff.com/suredry-scope-cabinets/>

Toll Free: (800) 433-4064 • [medicalsales@harloff.com](mailto:medicalsales@harloff.com)  
[www.harloff.com](http://www.harloff.com) • Ph: (719) 637-0300

**Optional dimmable LED lighting also available!**

# Why is a positive pressure system important?

The unique design of the SureDry Scope Cabinet includes both top fans for intake and lower fans for exhaust. This push/pull design creates a vortex of air that flows through the cabinet ensuring maximum drying efficiency.

## Intake:

Intake airflow comes through separate HEPA filters (which remove 99.7% of particulates, .03 Microns) driven by fully redundant fans producing a positive pressure vortex from the cabinet top. This creates an even circular motion of airflow over all scopes.

## Exhaust:

The lower exhaust fans pull the circulated air out of the cabinet ensuring the flow of air is consistent from top to bottom, providing a complete drying system for the entire scope. This push/pull multiple fan system (an industry first) ensures a positive pressure system to meet all major guidelines and it ensures forced drying and equal drying levels across the full length of each and every scope.



## SureDry Cabinets: A Smart Choice

to Meet and Exceed Current AORN, SGNA and AAMI Guidelines

Scope Storage Guidelines	Forced Air Drying	HEPA Filtered Air	Hang Scopes Vertically	Scopes Must Hang Freely, No Touching or Coiling	Easy to Clean and Disinfect	Closed Door Cabinet	Positive Air Pressure	Circulated Air Environment	Connectors — Use Disposables or Reusables Stored w/Scopes	Dry Scopes Completely	Uses Channel Drying (with Dri-Scope Aid®)
ANSI/AAMI ST91:2021		SureDry scope cabinets, with or without Dri-Scope Aid®, meet and exceed updated guidelines									
AORN	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
SGNA	✓	✓	✓	✓	✓	no data	✓	no data	✓	✓	✓
AAMI	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓



# 10 Reasons to Choose **HARLOFF** for Your Scope Cabinets



**Industry Best Positive Pressure HEPA Filtered Drying System.** Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

1

*A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines*

## Proudly Made in the USA



Harloff products are proudly **Designed, Engineered and Manufactured in the United States** at our Colorado Springs, CO facility.

6




## LED Lighting Option

**Brighten your cabinet** and more easily inspect scopes with fully dimmable optional overhead LED lighting which triggers on/off when you open and close the door.

2

## Easy Access to HEPA Filters

Change HEPA filters **without opening the cabinet.** Save time and ensure no cycle interruptions.

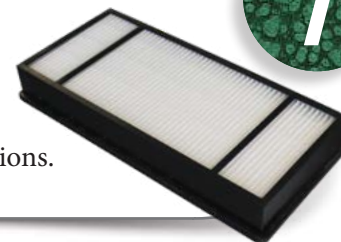

7

## Disposable Connectors

**Save time and ensure compliance!** Use single-use connectors for all internal channels of the endoscope. Part of the Dri-Scope Aid® package.



3

## Keyed or Keyless Locks

**Options include** key locks, basic electronic push button locks, advanced keypads and proximity readers. View locking options online at [www.harloff.com/locking-options/](http://www.harloff.com/locking-options/)



8

## Hooks and Holders for Specialized Scopes

Major brands to include **Fuji, Olympus, Pentax, Karl Storz and others.**



4

## Finish Options and Colors

Stainless Steel, Heat Sealed Powder Coat (High Gloss or Textured Low Gloss inside & out, view online at [www.harloff.com/finish-options/](http://www.harloff.com/finish-options/)) or upgrade to an **antimicrobial finish.**

9

## 12 Year Warranty

**Our Manufacturer's Warranty is Unmatched** in the industry. We build our cabinets to last and we stand behind them.

5

## Mobile Scope Cabinet Kit

Convert any size of scope cabinet into a mobile version. Includes 4 3" casters and steel mounting plate.



10