

FOR BRONCHOSCOPES

ENT scopes & other small scopes



ureDry
SCOPE CABINETS

**FOUR BRONCHOSCOPE OR EIGHT
ENT SCOPE CAPACITY WALL MOUNT
DRYING CABINET, KEY LOCK**

Scopes not included

SCW2448DRDP


This new wall mounted bronchoscope cabinet provides the same safe, efficient drying capabilities as that of scopes in other parts of your facility

- HEPA filtered Positive Pressure with two High Volume intake fans and one exhaust fan
- Customizable scope holders to fit most common scope types
- Holds up to 4 Bronchoscopes or up to 8 ENT Scopes
- Glass paneled doors for easy viewing of scope inventory
- Adjustable timer (10 to 180 mins) for fan operation.
- Fans activated automatically upon closing of doors
- Industry best 12 year limited manufacturer warranty
- Cabinet dimensions: 48" H x 24" W x 12.75" D



**Engineered &
Manufactured
in the U.S.**

**This optional LED-UV disinfection strip
helps to keep scopes sterilized!**

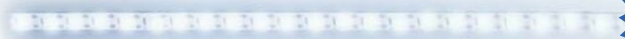

UV


HARLOFF
INNOVATING STORAGE FOR CARE

UV-C disinfection range is from 185 to 280nm. Our LED runs at 270nm, close to the optimal disinfection wavelength.

View similar products online at

<https://www.harloff.com/product-category/ent-and-bronchoscope-cabinets/>

**Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300**

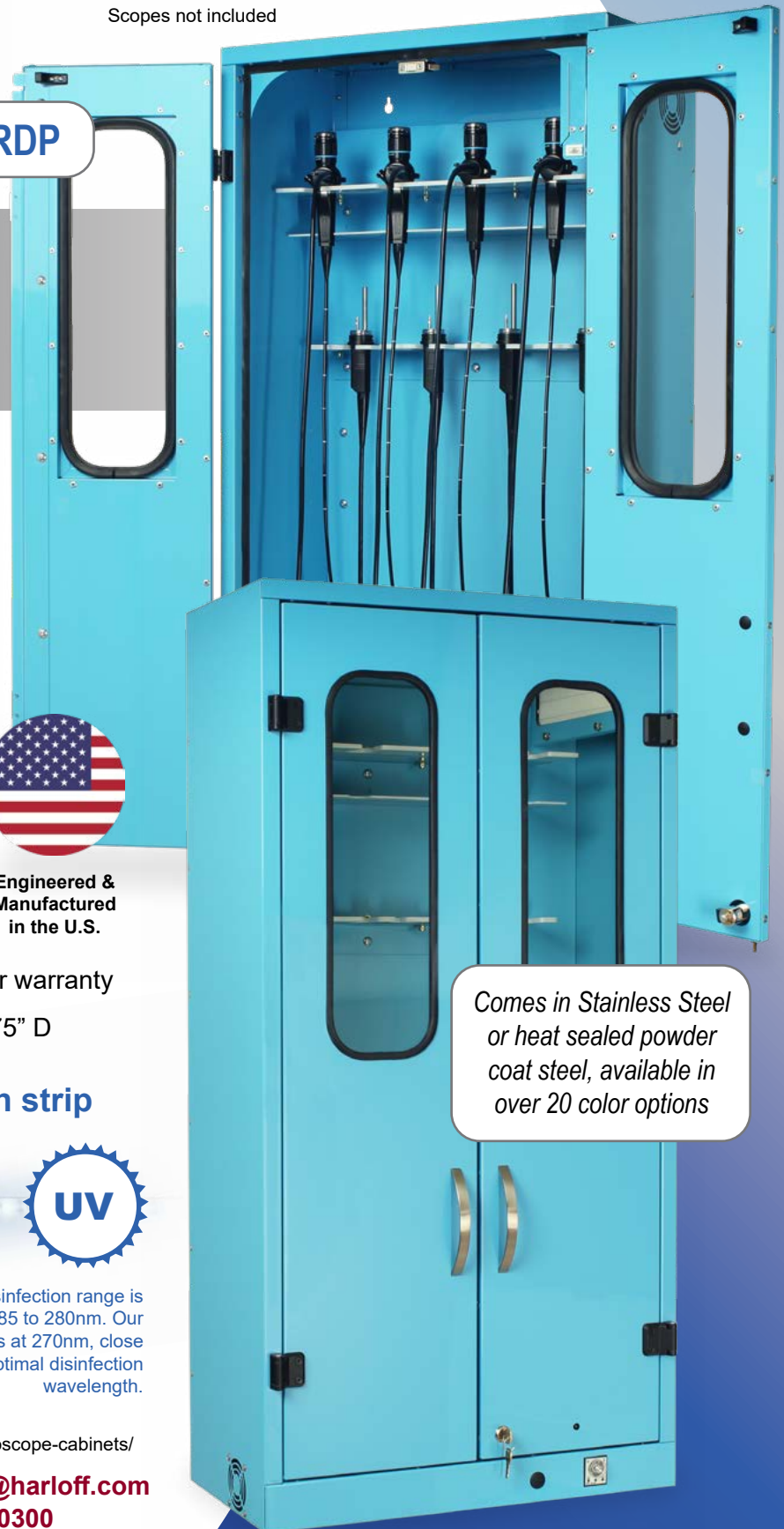

*Comes in Stainless Steel
or heat sealed powder
coat steel, available in
over 20 color options*

10 Reasons to Choose **HARLOFF** for Your Scope Cabinets

A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines



ureDry
SCOPE CABINETS

1 *Industry Best Positive Pressure HEPA Filtered Drying System.* Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

6 **Proudly Made in the USA**




Harloff products are proudly *Designed, Engineered and Manufactured in the United States* at our Colorado Springs, CO facility.

2 **Optional LED UV-Disinfecting Lights**

Help to keep your scopes sterilized with an optional UV radiation light strip as a way to disinfect your surfaces. Lights trigger in tandem with the fans after the cabinet doors are closed for scope drying.

7 **Easy Access to HEPA Filters**

Change HEPA filters easily with just thumb screws. Save time and minimize cycle interruptions.

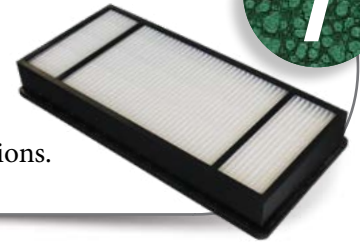

3 **Optional Wall or Cart Mounting**

Our cabinets fit your floor plan! Use our wall-mounted cabinet for permanent installations or choose a mobile cart-mounted cabinet to move scopes from your reprocessing area directly to your treatment room.

8 **Keyed or Keyless Locks**

Options include latch locks, key locks or standard electronic push button locks. View your locking options online at www.harloff.com/locking-options/



4 **Versatile Scope Holders**

Holds major brands including *Optim, Fuji, Olympus, Pentax, Karl Storz and others.*



9 **Scope Cabinet Finish Options**

Choose Stainless Steel, Heat Sealed Powder Coat in High Gloss or Textured Low Gloss, or upgrade to an *antimicrobial finish.*

5 **12 Year Warranty**

Our Manufacturer's Warranty is Unmatched in the industry. We build our cabinets to last and we stand behind them.



10 **Superior Color Choices**

Over 20 Colors to choose from. Cabinets have a heat sealed powder coat finish inside and out. View color options online at www.harloff.com/finish-options/